



Medical Refrigerator 2 to 20°C MR-1310



\$2,199.00

Laboratory Refrigerator	Laboratory Refrigerators and Freezers	Medical Refrigerator	Medical Refrigerator 2 to 20°C MR-1000
-------------------------	---------------------------------------	----------------------	--

Medical Refrigerator 2 to 20°C MR-1310 embedded with evaporator copper tube in a foam layer to prevent corrosion, to ensure heat transferring efficiency and glass door, different designed compact finned condenser for better efficiency. The maximum capacities of 317 L and defrost automatically. Build-in inner fan-driven and forced-draft air circulation system provides excellent temperature uniformity. Build-in microprocessor controller to certify tolerances of internal temperatures. Equipped with buzzer and flashing lights, power failure alarm, doors open alarm, high and low-temperature alarm.

Add To Cart


Product Description

Medical Refrigerator 2 to 20°C MR-1310 embedded with evaporator copper tube in a foam layer to prevent corrosion, to ensure heat transferring efficiency and glass door, different designed compact finned condenser for better efficiency. The maximum capacities of 317 L and defrost automatically. Build-in inner fan-driven and forced-draft air circulation system provides excellent temperature uniformity. Build-in microprocessor controller to certify tolerances of internal temperatures. Equipped with buzzer and flashing lights, power failure alarm, doors open alarm, high and low-temperature alarm.

Similar Products



Buy Now

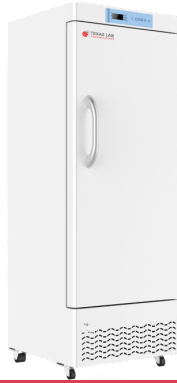

Buy Now



Buy Now



Buy Now